


MaxiGuard™ PRS

Patient Reminder System

Many times high-risk patients can become disoriented and endanger themselves by coming too close to the sides of their mattress. Today, in many long term care facilities, bed rails can no longer be used, which only compounds this danger.

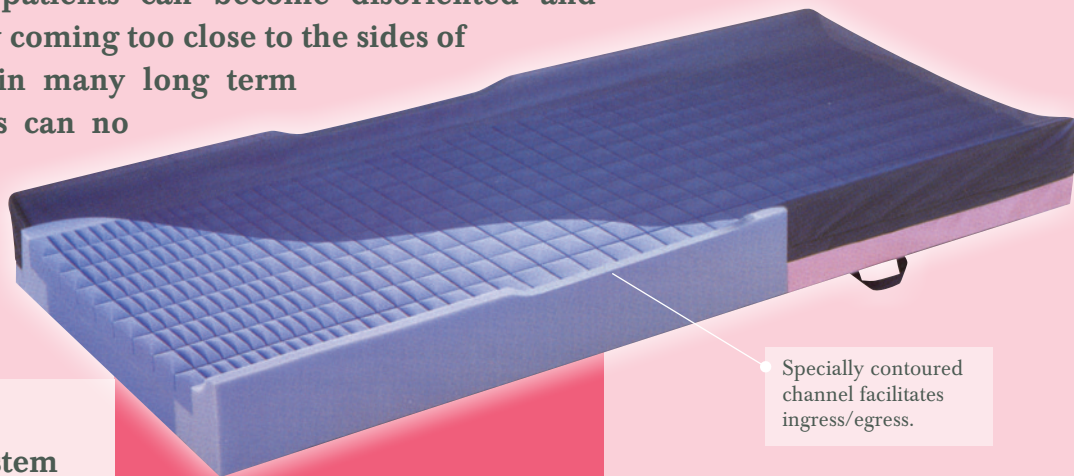


The New MaxiGuard Patient Reminder System

...incorporates raised *and* contoured side bolsters to reduce the risk of patients inadvertently falling out of bed. And, the MaxiGuard PRS is the only system which employs the time-tested and clinically proven floatation technology exclusive to MaxiFloat® and BG Industries; two names that you have trusted for decades for pressure management.

Other products currently on the market consist of little more than gluing angled borders onto large overlays. They lack the stability and resilience that only the PRS can provide.

Raised & contoured high density side bolsters

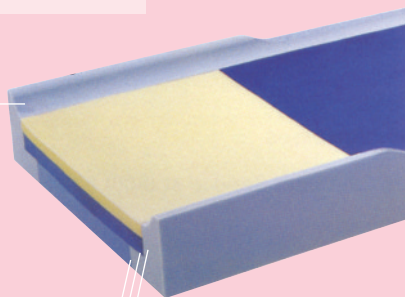

Boundary is bonded to side and bottom surfaces

MaxiGuard PRS option shown here on the MaxiFloat L series

The Unique Perimeter Boundary

By employing a high density foam boundary that is bonded not only to the sides of the mattress but to the bottom surface, the PRS helps keep your patient centered on the therapeutic laying surface.



Even in acute care facilities, the PRS can serve as a gentle reminder to your patients when they may be nearing the hard and unforgiving side rails attached to the bed frame.

Now both patients and caregivers can "rest" assured knowing that potentially dangerous falls or compromised patient positions between the mattress and bed rails can be minimised.

The BG Exclusives

Only the MaxiGuard PRS enlists all the features that you have grown used to with BG Industries:

- Hypolex® Top Cover with anti-microbial Vinyzene® protection on both sides.
- The patented anti-microbial flap behind the zipper.
- Our patented Linen Lock™ System holds bed linens in place.
- Heavy duty 360° YKK Zipper.
- Proven clinical outcome studies from both acute and long term care facilities.

The MaxiGuard PRS

The new patient reminder system from the people you have trusted for over three decades for therapeutic pressure management.

The MaxiGuard Patient Reminder System

Safe and Gentle Protection for Your High Risk Patients.



SEXTON

Sexton Trading Company Pty. Ltd.

ABN 68 009 829 485

72 Crockford Street, Northgate QLD 4013 Telephone: (07) 3265 9999
PO Box 135 Hamilton Qld 4007 Australia Web: www.sexton.com.au

Facsimile: (07) 3260 6079
Email: sexton@sexton.com.au




Quality Management System
Certificate Registration No. 1058
AS/NZS ISO 9001: 2008

MaxiGuard™ PRS

Patient Reminder System

Specifications:

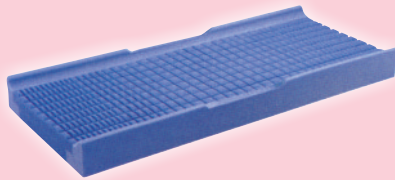

The Inner Foam Core:

The top layer is medium density foam, contour cut in seven zones along the full length of the mattress. The protective foam side rails of the PRS rise 2 inches above the top layer.

The bottom support layer consists of high density foam, 28" wide and 3" thick. 3½" wide contoured foam side rails are bonded to the edges.

All the components of the inner foam core meet or exceed the specifications in California Upholstery Provisions as stated in Technical Bulletin #117.



The Top Cover:

The Hypolex® top cover offers very low moisture vapor transmission rates and is anti-microbial. The top cover can be cleaned with standard cleaning procedures, using properly diluted solutions. The Hypolex cover system can be removed and replaced, if necessary and is attached with a 360° YKK Zipper.

The Bottom Cover:

The bottom cover is an 11 oz., 1,000 denier aberdeen rose vinyl that is both waterproof and anti-bacterial. It includes our patented anti-microbial flap behind the zipper preventing the ingress of fluids or other foreign material directly into the mattress interior. Four reinforced turning handles are included.

Ordering Information:

Name	Size	Order Number
MaxiGuard PRS Mattress	35 x 75 x 6 35 x 78 x 6	

Warranty:

This warranty applies only to defects in structure, workmanship and materials and only to damages arising from normal hospital usage on a non pro rated basis for a period of three (3) years from date of purchase. Damages arising from abnormal use such as those caused by needle punctures, burns, chemicals, negligent use or improper care or cleaning, are excluded from coverage. Manufacturer will pay shipping and handling charges incurred in connection with this warranty during the first year; thereafter, those charges are the facilities' sole responsibility.

The MaxiGuard PRS is also available on our MaxiFloat mattresses.
If you have any questions or require further information, please contact:



Sexton Trading Company Pty. Ltd.

ABN 68 009 829 485

72 Crockford Street, Northgate QLD 4013 Telephone: (07) 3265 9999 Facsimile: (07) 3260 6079
PO Box 135 Hamilton Qld 4007 Australia Web: www.sexton.com.au Email: sexton@sexton.com.au



Quality Management System
Certificate Registration No. 1058
AS/NZS ISO 9001: 2008